



Narrow Head Rest


Description

A T-shaped head rest which can be used in place of a standard OR table head rest, granting the surgeon more easily access to the patient's head during ENT, neurosurgical, ophthalmic, and other such surgeries. Its radiolucent construction enables it to be used for patient imaging, and its telescoping mounting rails allow it to fit any OR table with the use of two blade-style sockets (sold separately; see *PRODUCT COMPATIBILITY*).

Reusable pressure management pad included.

Made in USA.

MPR-800-0052

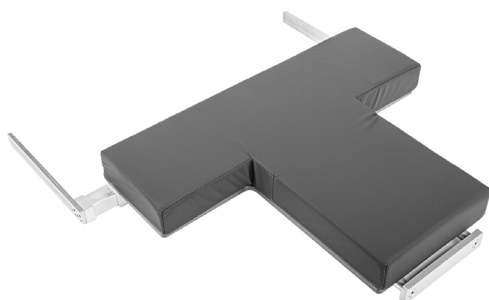

Product Specifications

Head Rest Length	16" (40.64cm)
Mounting Column Length	9.25" (23.5cm)
Base Width	20" (50.8cm)
Head Width	8" (20.32cm)
Pad Thickness	2" (5.08cm)

Product Compatibility

#	DESCRIPTION
MPR-BC	MPR Basic Blade Side Rail Socket
MPR-SRS	MPR Aluminum Side Rail Socket
MPR-800-0228	Simple Clamp

For Further Information

If further information on this product is required, please contact MPR Orthopedics Customer Service at 800.524.5194 in the United States, or your authorized representative, or visit the MPR Orthopedics website at www.mprortho.com.