MPR
ORTHOPEDICS

MPR Arthroscopic Stress Post

MPR-9305



Instructions for Use

MPR Arthroscopic Stress Post (MPR-9305)

Description

The MPR Arthroscopic Stress Post is a positioning device primarily used during knee arthroscopic procedures to provide leverage of the leg, granting better access to the knee medial compartment.

The device requires no assembly before use, and consists of a waterproof durable molded polyurethane pad and stainless steel mounting post, the height of which can be freely adjusted to meet the needs of the procedure. The device can also be draped during procedures, negating the need for any additional covering or sterilization and making it easy to clean and disinfect.

Compatible with any ⅝" Clark socket, such as the MPR-CS, MPR-SRC, or MPR-CS-AIO (*sold separately; see ACCESSORIES*).


Indications

For use during knee arthroscopic procedures to improve access to the knee medial compartment.

Cautions

- UNITED STATES FEDERAL LAW RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A PHYSICIAN.
- This product is provided non-sterile and must be properly cleaned prior to use. Refer to the CLEANING AND MAINTENANCE section on page 2 of this document.
- Prior to use, examine the instrument for possible damage and to ensure proper functioning.
- A surgeon *should not* attempt clinical use of this device before reviewing these instructions for use.

Components

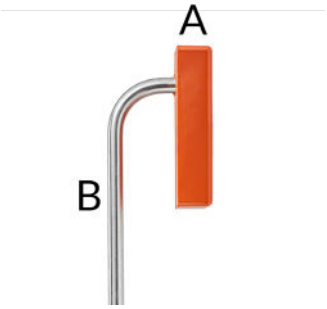

A	Face Plate w/ Polyurethane Pad
B	Mounting Post

Accessories

REF	Description
MPR-CS	MPR Stainless Steel Clark Socket
MPR-SRC	MPR Over Side Rail Clark Socket
MPR-CS-AIO	MPR All-In-One Clark Socket
MPR-URC	MPR Aluminum Universal Rail Clamp
MPR-ORUC	MPR Over the Rail Universal Clamp

MPR ORTHOPEDICS

MPR Arthroscopic Stress Post (MPR-9305)

Instructions for Use

1. Position patient on operating table so that they are lying supine with their feet at the end of the bed and ipsilateral hip at the edge. Apply a tourniquet to the patient's upper thigh.

2. Attach a ⅝" Clark socket or equivalent to the table rail.

3. Insert the mounting post into the Clark socket and secure. Keep it as upright as possible.

4. Position the stress post around 5cm superior to the proximal pole of the patient's patella. It should be positioned so that when the patient's thigh is pressed against the stress post, you have enough space to stand between the patient's ankle and the operating table.

TO EXPOSE THE MEDIAL COMPARTMENT

5. Have either you or your assistant grasp the patient's foot and push it laterally into valgus, holding it so that the thigh is pressed against the stress post. If you do not have an assistant, you can brace the leg against the stress post with your hip while operating.

Cleaning & Maintenance

1. Preventing blood, saline, or other materials from drying on the device will make cleaning easier. *Do not* immerse the device for cleaning.

2. Manually clean the device with a soft bristle brush and enzyme solution prepared according to the manufacturer's recommendations. Rinse with tap water.

3. Apply a surface disinfectant spray to the device. Let stand for the time recommended by the manufacturer.

4. Rinse the device thoroughly in lukewarm, running tap water for one minute.

5. Place the cleaned device on a sheet of absorbent matting and allow it to air dry.

⚠ ***Do not* use an autoclave to sterilize the device.**

Warranty

MPR Orthopedics products are guaranteed to be free from defects in material and workmanship for the warranty period for a particular product, *beginning from date of invoice*. Contact MPR Orthopedics for specific warranty information.

This limited warranty is restricted to repair or replacement by MPR Orthopedics, at its option, of any product found to be defective during the warranty period. Damage inflicted on a product by the user that causes it to be unsuitable for refurbishment may result in additional charges, regardless of warranty status. All warranties apply to the original buyer *only*.

In no event shall MPR Orthopedics be liable for any anticipated profits, consequential damages, or loss of time incurred by the buyer with the purchase or use of any product.

NO OTHER WARRANTY, EXPRESSED OR IMPLIED, IS GIVEN.

****WARRANTY IS VALID FOR 1 YEAR***

MPR Orthopedics Offers:

- Free trials on stress posts.
- Repairs for defective stress posts.
- Trade-ins for excess equipment.
- Clark sockets.
- Universal clamps.

For Further Information

If further information on this product is required, please contact MPR Orthopedics Customer Service at 800.524.5194 in the United States, or your authorized representative, or visit the MPR Orthopedics website at www.mprortho.com.

Also Available

REF	Description
MPR-5600	MPR Stress Post
MPR RT.0010	MPR Shutt Universal Lateral Stress Post
MPR RT.0009	MPR Shutt XL Universal Later Stress Post
MPR RT.0010NV	MPR Shutt Universal Lateral Stress Post w/o Velcro
MPR RT.0009NV	MPR Shutt XL Universal Later Stress Post w/o Velcro
MPR-8840	MPR Leg Stabilizer

A stable, time-proven orthopedic product.